



OPHTEISBIO 1.6

Cohesive

Biofermented NaHA 1.6%

OPHTEISBIO 1.8

Cohesive

Biofermented NaHA 1.8%

OPHTEISBIO 3.0

Dispersive

Biofermented NaHA 3.0%



All of our OphteisBio products

- have optical clarity for maximum visibility.
- are made with a Biofermented Sodium Hyaluronate.

OPHTEISBIO 1.6%

Cohesive • Designed for all types of surgery

- Optimal maintenance of volume in the anterior chamber or capsular bag.
- Good coating power.
- Easy injection and removal during surgical stages.

OPHTEISBIO 1.8%

Cohesive • Designed for small incision


- Good cohesivity at low shear-rate for a stable anterior chamber.
- Strong coating of tissue thanks to improved dispersive property.
- Easy to aspirate with high molecular weight.

OPHTEISBIO 3.0%

Dispersive • Designed for excellent endothelial protection

- Low molecular weight, high NaHA concentration.
- Assures maximum protection and viscosity.
- Good maintenance of anterior chamber.

Rheological Profiles

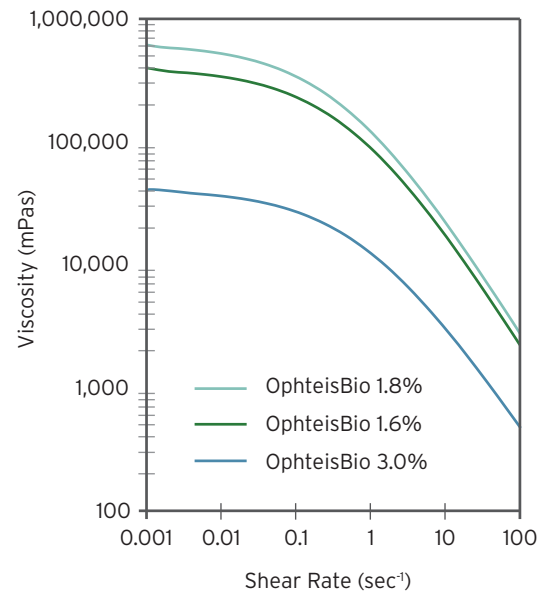

Product	OphteisBio 1.6%	OphteisBio 1.8%	OphteisBio 3.0%
Polymer Origin	Biofermentation	Biofermentation	Biofermentation
Sodium Hyaluronate Concentration	1.6%	1.8%	3.0%
Molecular Weight (Dalton)	approx. 3 million	approx. 3 million	approx. 0.75 million
Zero Shear Viscosity (mPas)	avg. 400,000	avg. 600,000	avg. 30,000
Osmolality (mOsm/kg)	300 to 350	300 to 350	300 to 350
pH	6.8 - 7.6	6.8 - 7.6	6.8 - 7.6
Shelf Life (years)	3	3	3
Storage	2 to 25°C	2 to 25°C	2 to 25°C
Syringe Volume (ml)	1.1	1.1	1.1
Cannula Gauge (G)	27	27	25



Regulatory restrictions may apply in certain markets.
©Rayner 2018. OphteisBio is a registered mark of Rayner.
EC 2018-02/12/18

GLOBAL HQ & UK OFFICE
RAYNER INTRAOCULAR LENSES LIMITED
10 Dominion Way
Worthing, West Sussex
BN14 8AQ, United Kingdom

Tel: +44 (0) 1903 258 900
Email: iol_enquiries@rayner.com
Web: rayner.com

UK CUSTOMER SERVICE
Tel: +44 (0) 1903 854 560
Email: orders@rayner.com

GERMANY CUSTOMER SERVICE
Tel: +49 30 6290783 - 0
Email: info.de@rayner.com

ITALY CUSTOMER SERVICE
Tel: +39 800 972 485
Email: ordini@rayner.com

SPAIN CUSTOMER SERVICE
Tel: +349 1048 0055
Email: pedidos@rayner.com

PORTUGAL CUSTOMER SERVICE
Tel: +351 221 206 400
Email: pedidos@rayner.com